


C-PAP


What Is CPAP?

CPAP (Continuous Positive Airway Pressure) is a device that lets you breathe better at night. CPAP can help if you have sleep apnea. This serious health problem can cause you to snore, gasp, or snort in your sleep. Snoring can keep you and your partner awake at night. You may also feel tired in the morning. Talk to your health care provider to find out if CPAP is right for you.

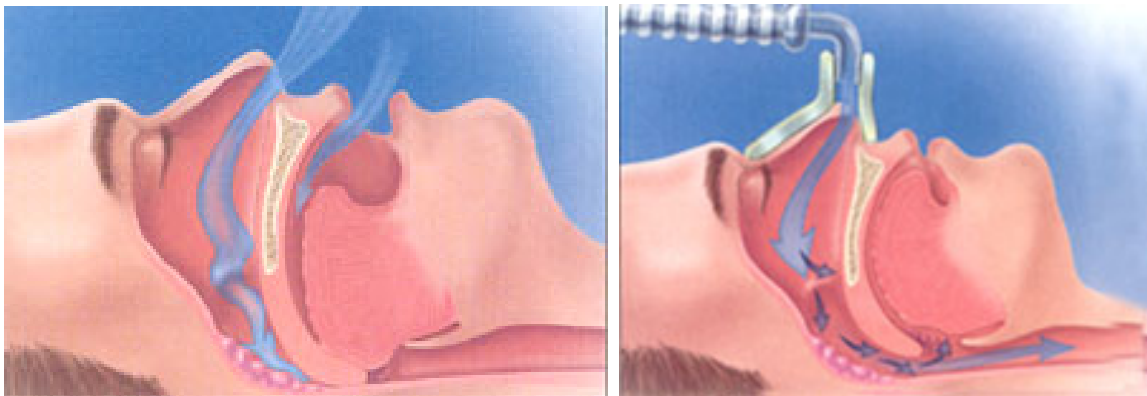

What Is Sleep Apnea?

People with sleep apnea stop breathing for at least 10 seconds at a time while sleeping. With a loud gasp, breathing begins again. This pattern is repeated all night, often with sudden body movements and irregular heartbeat. Sleep apnea can cause morning headaches, mood changes, and memory problems. It can also be a risk factor for heart disease or stroke. Sleep apnea may get worse with age or weight gain.

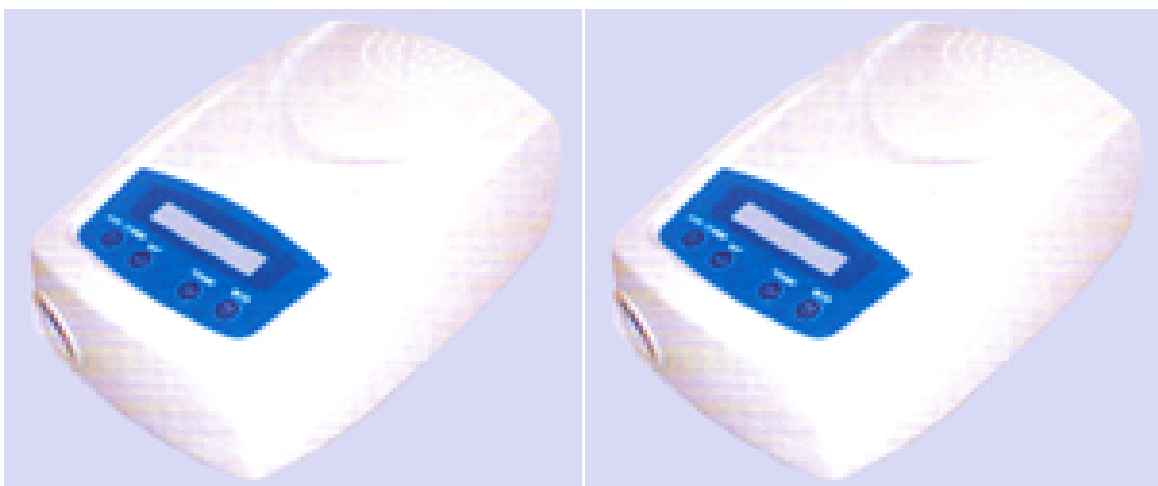

CPAP Can Help

Your health care provider can determine if you have sleep apnea. Then he or she may recommend CPAP. The CPAP device has a soft plastic mask that fits over your nose. The

mask is connected to the unit by a flexible tube. When you go to bed each night, you place the mask over your nose and turn on the machine. The machine sends a gentle, steady stream through your airway to keep it open. This allows you to breathe freely.

Features:

Easy to-carry design

It is compact and easy to carry with small size and integrated hidden handle on the backside. Traveling with CPAP is no longer troublesome.

LCD display with back light and easy-to-use interface

The LCD display is easy to read, and the intuitive user interface is simple to use.

Pressure level can be adjusted based on prescription

Operating pressure, ranging from 4 to 18cm H₂O, can be pre-set to company with physician and closely matched with individual needs.

Ramp time can be set for delay delivering pressure

Deluxe ramp capability allows a 0-45 minute adjustable time period to help user fall to sleep easier.

Compliance meter for reimbursement support

The therapy hours can be displayed in details at patient usage for reimbursement to therapeutic purposes

Extremely quiet operation

Specification:	
Pressure Range	4 - 18 cmH ₂ O (1 cmH ₂ O increment)
Default Pressure	4 cm H ₂ O
Ramp Time	0 - 45 minutes (5 min increment)
Sound Output	<30dBA
Control	START/STANDBY, RAMP, PRESSURE, METER
Indicator	1 x 8 character LCD

Electrical Characteristics	Universal power supply 85-255 VAC, 50/60 Hz (47-63 Hz) 110/230 VAC Power Cord Option
Size (WxDxH)	17.1cm x 25.3cm x 11cm
Operating Temperature	+ 5oC to + 40oC/ + 41Fo to +104Fo
Environment Humidity	15%-95% non-condensing